**MiniSpacer\* Dual Spray Metered Dose Inhaler Adapter**

Bi-Directional Spray and Cap

15 mm O.D. / 15 mm I.D. ISO Connections

Accepts plastic and metal pMDI canister stems, including canisters with integrated dose counters.

#### INDICATIONS FOR USE

Model number 1543 is a single patient, disposable device for dispensing pressurized metered dose inhaler (pMDI) medication into a breathing circuit, as prescribed by a physician or other licensed health care practitioner. The device is indicated for patients on a breathing circuit, for whom aerosol medication has been prescribed, in short and long term critical care environments. REF 1543 is intended for use only when connected to a 15 mm fitting.

#### DIRECTIONS FOR USE

**Installing the MINISPACER Adapter in the Patient Circuit**

- Remove the MINISPACER adapter from package.
- Install the MINISPACER device in the patient circuit as indicated here or according to facility guidelines:
  - For a 22 mm circuit, install the MINISPACER between the patient wye and endotracheal tube.
  - For a 15 mm circuit, install the MINISPACER between the patient wye and the inspiratory limb of the patient circuit.

**CAUTION** : Ensure the device is securely installed before use.

**NOTE** : When correctly installed, the MINISPACER MDI port will be vertical and the pMDI canister will point downward during treatment.

- Place the cap firmly on the MINISPACER port.

**Administering MDI medication via the MINISPACER Adapter**

- Uncap the MINISPACER MDI port.
- Insert the MDI canister stem into the MINISPACER port.
- Coordinate MDI actuation with the beginning of inspiration.

**CAUTION : Support the MINISPACER adapter while depressing the canister.**

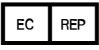
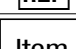
- Administer the dose as ordered.
- After medication has been delivered as prescribed, remove the pMDI canister and firmly cap the MINISPACER port.
- Monitor the patient response.

**NOTE** : Discard the MINISPACER adapter when replacing the circuit according to facility guidelines, not to exceed two weeks.

**CAUTION:** REF 1543 MINISPACER adapter will not increment the counter on pMDI canisters with integrated dose counters surrounding the stem.

**WARNING** : Do not clean and reuse or sterilize. Potential risks associated with cleaning and reusing, reprocessing and/or sterilization include: cross-contamination or infection resulting from inadequate cleaning, product degradation, and/or failure of the device to perform as intended.

#### Symbol Glossary

	Manufacturer		Consult instructions for use
	Date of manufacture		Authorized representative in the European Community
	Batch code/lot number		Prescription only. U.S. Federal law restricts this device to sale by or on the order of a physician or properly licensed practitioner.
	CE Mark		Caution. The symbol can also be used to mean Attention, see Instructions for Use.
	Catalog, reorder or reference number		

Item No.	Description
 1543	MiniSpacer® Dual Spray MDI Adapter, 15 mm O.D. / I.D.

	Manufacturer		0482
	<b>Thayer Medical Corporation</b> 4575 S. Palo Verde Rd. Tucson, AZ 85714-1961 USA www.thayermedical.com Made in USA		MDSS GmbH Schiffgraben 41 30175 Hannover, Germany

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The MINISPACER dual-spray adapter was evaluated for particle size and dosing characteristics using an 8 stage cascade impactor. The MINISPACER was connected to the cascade impactor with ventilator tubing, wye and an endotracheal (ET) tube as they would be used in a ventilator circuit. Air was drawn through the assembly at a flow rate of 28.3 l/min. The aerosol samples exited the ET tube into the cascade impactor and comprised the Total Delivered Dose.

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